




National and International Standards demand that vacuum-assisted steam sterilizers must be tested at the beginning of each working day. Test Pack designed for vacuum-assisted sterilizers to test the efficiency of the vacuum system, to remove air from the sterilizer chamber and to detect leaks in the system that would allow air into the chamber. A uniform change of the indicator sheet indicates that all the air was removed and replaced by steam it also indicates that steam penetration occurs in a healthy way. Performed after a warm-up cycle at the beginning of every day, the Bowie & Dick Test verifies the sterilizer effectively removes air.

- ★ It is easy to use and cost effective.
- ★ Water-based and non-toxic chemical indicator is used.
- ★ Indicator pattern is easy to read and interpret.
- ★ Designed to give equivalent results to those obtainable by the Bowie & Dick protocol.
- ★ The backside of the indicator sheet contains fields to be filled out and it can be kept as a permanent record for traceability purposes.
- ★ Compatible with ISO 11140-4 (Type 2 + 4) (Class 2 + 4).
- ★ Different designs are made for the company/person on request.

BD-200

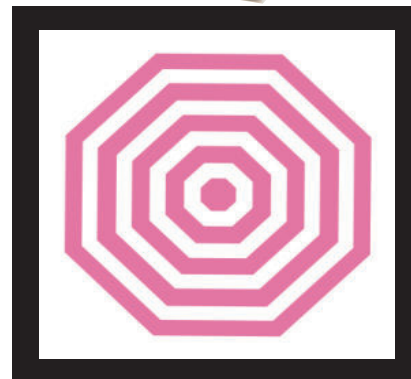

CARD SIZE:
124 X 124 mm

REF NO	DESCRIPTION	PACKAGING
BD-200	Bowie & Dick Test Pack (200 Sheets)	80 Pcs. / CARTON
BD-400	Bowie & Dick Test Pack (400 Sheets)	48 Pcs. / CARTON