


The Ideal Surgical Solution for
Ophthalmology, Plastic Surgery, ENT, and
Outpatient/Same Day Surgery

Hausted®
SurgiChair Series



HAUSTED®
PATIENT HANDLING SYSTEMS, LLC

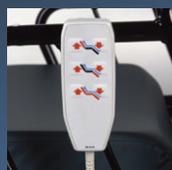

*Model ESC-EYE
with electric/battery
powered controls.*

OPERATIONAL BENEFITS

- The articulating head piece with dual operating control knobs allows infinite adjustability and precise microscopic movement of the head section.
- Low profile head piece design offers excellent site access for both superior (over brow) and lateral (side approach) procedures.
- The chair offers simple controls for adjustment of height, backrest, leg section, and Trendelenburg.
- A wide variety of surgical accessories can be adapted to the unit to provide versatility.

ECONOMIC BENEFITS

- Pre-induction, transport, treatment, and recovery can be performed on the same unit to improve case turnaround and streamline patient handling.
- Reduce staff time and prevent staff injury involved in transferring patient during various phases of care.



Pre-Induction/Transport



Converts easily to surgical position



Post-Anesthesia Care/Recovery



PATIENT OUTCOME AND SAFETY BENEFITS

- With the hand held control, patients can position themselves during recovery, which minimizes staff involvement and can help shorten recovery time.
- A full range of patient positioning is achieved to allow patient comfort during every phase of care.



Microscopic articulation of the head piece is achieved by rotating the control knobs to raise/lower the patient's forehead or chin. The patient's head is continually supported during head articulation to provide optimal safety during all positioning.

STANDARD FEATURES

- 5" (12.7 cm) easy rolling casters
- Four wheel brake and steer system operable from both sides of unit
- Easy to clean thermoformed base cover
- Polymer coated steel litter top
- Low profile articulating headpiece
 - *Dual operating control knobs*
- Fold-down foot rest (height adjustable)
- Pre-Op/Post-Op Head Extension
- Accessory wells (2 on each side)
- Padded armrests maintain patient security
- 3" (7.6 cm) comfort mattress and head pad


STANDARD FEATURES FOR ELECTRIC/BATTERY POWERED 5E8-EYE-ST

- Power controls for:
 - *Height adjustment*
 - *Backrest*
 - *Leg section*
 - *Trendelenburg*
- Patient/Staff Hand Held Control
- Emergency Backrest Release lever to override power backrest
- Battery back-up allows power functions when not plugged in

OPTIONAL ACCESSORIES

- Telescoping IV Pole
- Oxygen Tank Holder
- 4" (10 cm) PressureCare mattress
- Concave support in 2" (5 cm), 3" (7.6 cm), or 4" (10 cm) thickness


- Adjustable wrist rests
 - *Dual side (lateral/temporal approach)*
 - *Full "U" (superior/over-brow approach)*
- Universal patient/instrument tray
- Fold down side tables
- Armboard with 2" (5 cm) pad and clamp
- Surgical accessory rail(s)
- Oxy-Flex oxygen delivery system/drape support
- Foot control pedal for electric/battery powered hi/lo
- Knee crutch/leg holder with mounting adapters
- Heel stirrups with mounting adapters
- Paper roll holder
- Patient security straps with velcro or buckle
- Five custom color pad options: Gray, Tan, Mauve, Blue, Green, Black (Standard)

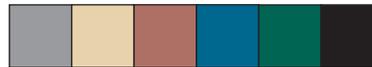

GENERAL SPECIFICATIONS

MODELS	ESC-EYE-ST ESD-EYE-ST*
Height Range: High	31.5 in [80 cm]
Height Range: Low	23.5 in [60 cm]
Overall Width	28 in [71 cm]
Overall Patient Surface Width	22 in [56 cm]
Overall Length	75 in [191 cm]
Overall Patient Surface Length	75 in [191 cm]
Backrest Degree of Movement	0-80°
Standard Trendelenburg	6°
Siderail Height	12 in [30.5 cm]
Siderail Length	29 in [74 cm]
Maximum Weight Capacity	325 lbs [147.5 kg]

*230 VAC (+/-10%)

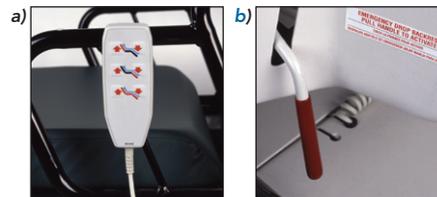

Electronic Controls

Safety and convenience are offered with the Electronic controlled version including (a) hand held control, (b) emergency backrest release lever



PATIENT HANDLING SYSTEMS, LLC

334.215.5151 **Main**
 877.706.5151 **Toll Free**
 334.215.5150 **Fax**

2511 Midpark Drive
 Montgomery, AL 36109

www.hausted.com