

## STERILISATION PACKAGING

Our flexible production methods ensure that the best sizes of wrap can be produced for our Customers, even if these are not industry standard sizes, without compromising performance or protection levels.

Our Barrier Wraps provide a bacterial barrier and protect the contents of a set or pack from contamination. When classed as alcohol repellent, barrier wraps can be used within the sterile field.

Choice Wrap is designed to reduce the risk of rips, tears, and wet packs. Choice SMMS are made of polypropylene and come in three weights: 43gsm, 55gsm, 70gsm. Any two of the Choice Wraps can be interleaved in blue or green to meet the sterilisation service user requirements.

Want to trust your Transport Wrap? All Transport Wraps must be used in conjunction with one compliant Barrier Wrap to completely prevent the ingress of micro-organisms. Use of a single barrier layer with Transport Wrap can be a cost-effective solution that reduces damage in handling and transport. If adopted, changes to both wrapping and unwrapping procedures may be required to ensure aseptic opening.

Sterilisation Pouches and Reels are also available to ensure that the wrapping process is efficient and effective. Made to withstand vigorous processing and handling, pouches and reels are made from medical grade film with breathable webs and suit all common sterilisation methods.

Barrier Wrap

Crepe Wrap

Liners

Sterilisation Tape

Transport Wrap

Amcors Pouches and Reels

Wipak Pouches and Reels

Please note: Unless specified, all products are not made with natural rubber latex.





## Sterilisation Packaging – Barrier Wrap

### CLINIPAK CHOICE WRAP

A fully synthetic SMMS product designed to offer flexibility and enable users to pack conventionally or with 'two layers' at once. Clinipak Choice Wrap is strong and available with alcohol repellence in both green and blue for use within the sterile field. Wraps are also available with blue and green interleaved. Compatible with steam and EtO sterilisation. Single use.

#### A CLINIPAK CHOICE WRAP BLUE

Product code	Size	GSM	Qty
CB060	60 x 60cm	55gsm	350 sheets
CB075	75 x 75cm	55gsm	300 sheets
CB090	90 x 90cm	55gsm	300 sheets
CB100	100 x 100cm	55gsm	200 sheets
CB120	120 x 120cm	55gsm	120 sheets
CB1315	130 x 150cm	55gsm	100 sheets
CB1715	170 x 150cm	55gsm	50 sheets
CB1815	180 x 150cm	55gsm	50 sheets

#### B CLINIPAK CHOICE GREEN WRAP

Product code	Size	GSM	Qty
CG060	60 x 60cm	55gsm	350 sheets
CG075	75 x 75cm	55gsm	300 sheets
CG090	90 x 90cm	55gsm	300 sheets
CG100	100 x 100cm	55gsm	200 sheets

#### C CLINIPAK CHOICE BLUE LIGHT WRAP

Product code	Size	GSM	Qty
CBA050	50 x 50cm	43gsm	300 sheets
CBA060	60 x 60cm	43gsm	350 sheets
CBA090	90 x 90cm	43gsm	300 sheets
CBA100	100 x 100cm	43gsm	200 sheets
CBA120	120 x 120cm	43gsm	120 sheets
CBA1415	140 x 150cm	43gsm	100 sheets





## Sterilisation Packaging – Barrier Wrap

### A CLINIPAK CHOICE BLUE ULTRA WRAP

Clinipak Choice Blue Ultra is extra strong, drapeable, and available with alcohol repellence in both green and blue for use within the sterile field. Compatible with steam and EtO sterilisation. Single use.

Product code	Size	GSM	Qty
CBU1014	100 x 140cm	70gsm	120 sheets
CBU120	120 x 120cm	70gsm	120 sheets
CBU1315	130 x 150cm	70gsm	100 sheets
CBU140	140 x 140cm	70gsm	100 sheets

### B CLINIPAK CHOICE GREEN X WRAP

A fully synthetic SMMS product designed to offer flexibility and enable users to pack conventionally or with 'two layers' at once. Available with alcohol non repellent in green for use within the sterile field. Compatible with steam and EtO sterilisation. Single use.

Product code	Size	GSM	Qty
CGX090	90 x 90cm	55gsm	300 sheets
CGX100	100 x 100cm	55gsm	200 sheets

A



B



K



## Sterilisation Packaging – Barrier Wrap

### CLINIPAK CHOICE WRAP DRY

A patented tray wrap solution is designed to heavily reduce the occurrence of wet packs. The inner surface exhibits a high degree of absorbency and the outer surface repels both water and alcohol. The absorbent surface disperses the condensation within a typical steam sterilisation process and assists with evaporation during the vacuum cycle. This wrap is manufactured with a blend of materials that are compliant with the recognised standards for toxicology, bacterial resistance, and drape. Compatible with steam and EtO sterilisation. Single use.

### B CLINIPAK CHOICE WRAP DRY GREEN

Product code	Size	GSM	Qty
CDG075	75 x 75cm	79gsm	300 sheets
CDG090	90 x 90cm	79gsm	300 sheets
CDG100	100 x 100cm	79gsm	200 sheets
CDG120	120 x 120cm	79gsm	120 sheets
CDG140	140 x 140cm	79gsm	100 sheets
CDG150	150 x 150cm	79gsm	100 sheets
CDG1315	130 x 150cm	79gsm	100 sheets

B



K

Sterilisation Packaging



## Sterilisation Packaging – Barrier Wrap

### A RELIANCE® SOLO WRAP

A two-sheet wrap system, with two colours of latex-free wrap joined together along the sides. Each sheet is a five-layer nonwoven tray wrap, manufactured from spun-bonded and melt blown polypropylene, providing a wet and dry microbiological barrier. Bonding the wrap simplifies the wrapping process whilst increasing durability to avoid rips, tears or punctures. Wraps are also available interleaved with other compliant barrier layers, or transportation wrap. Blue and green bonded wrap. Compatible with steam and EtO sterilisation. Single use.

### RELIANCE® SOLO WRAP 200B/G

Product code	Size	GSM	Qty
CT43LGT-60	60 x 60cm	43gsm	125 sets
CT43LGT-75	75 x 75cm	43gsm	150 sets
CT43LGT-90	90 x 90cm	43gsm	50 sets
CT43LGT-100	100 x 100cm	43gsm	50 sets
CT43LGT-120	120 x 120cm	43gsm	50 sets
CT43LGT-140	140 x 140cm	43gsm	25 sets
CT43LGT-1418	140 x 180cm	43gsm	25 sets

### RELIANCE® SOLO WRAP 300B/G

Product code	Size	GSM	Qty
CT50MED-60	60 x 60cm	50gsm	125 sets
CT50MED-75	75 x 75cm	50gsm	50 sets
CT50MED-90	90 x 90cm	50gsm	50 sets
CT50MED-100	100 x 100cm	50gsm	50 sets
CT50MED-120	120 x 120cm	50gsm	50 sets
CT50MED-140	140 x 140cm	50gsm	25 sets
CT50MED-1418	140 x 180cm	50gsm	25 sets

### RELIANCE® SOLO WRAP 400B/G

Product code	Size	GSM	Qty
CT60HVV-60	60 x 60cm	60gsm	125 sets
CT60HVV-75	75 x 75cm	60gsm	50 sets
CT60HVV-90	90 x 90cm	60gsm	50 sets
CT60HVV-100	100 x 100cm	60gsm	50 sets
CT60HVV-120	120 x 120cm	60gsm	25 sets
CT60HVV-140	140 x 140cm	60gsm	25 sets
CT60HVV-1418	140 x 180cm	60gsm	25 sets





## Sterilisation Packaging – Crepe Wrap

### RELIANCE® NON-WOVEN WRAP

Non-Woven Wrap is a high-performance sterilisation material that is suitable for inner and outer wrapping of medium to large trays and packs. This is a wet laid wrap with 50% cellulose and 50% synthetic fibres that improve strength and increase softness. These products are also available in interleaved green and blue wrap. Single use.

	Product code	Colour	Size	GSM	Qty
A	RELBLU090	blue	90 x 90cm	60gsm	200 sheets
	RELBLU100	blue	100 x 100cm	60gsm	200 sheets
	RELBLU120	blue	120 x 120cm	60gsm	100 sheets
B	RELGRN100	green	100 x 100cm	60gsm	200 sheets
	RELGRN1315	green	130 x 150cm	60gsm	100 sheets

### RELIANCE® CREPE WRAP


Crepe paper is suitable for lighter instruments and sets. Manufactured from 100% medical quality cellulose fibres and combines most economical wrapping with good drapability. A soft crepe and easier to handle due to the creping process. Compatible with steam and EtO sterilisation. Single use.

	Product code	Colour	Size	GSM	Qty
C	GRNCREPE025	green	25 x 25cm	60gsm	2000 sheets
	GRNCREPE060	green	60 x 60cm	60gsm	500 sheets
	GRNCREPE075	green	75 x 75cm	60gsm	250 sheets
	GRNCREPE090	green	90 x 90cm	60gsm	250 sheets
	GRNCREPE100	green	100 x 100cm	60gsm	250 sheets
	GRNCREPE1014	green	100 x 140cm	60gsm	125 sheets
	GRNCREPE120	green	120 x 120cm	60gsm	125 sheets
D	WHCREPE1315	green	130 x 150cm	60gsm	125 sheets
	WHCREPE030	white	30 x 30cm	60gsm	2000 sheets
	WHCREPE050	white	50 x 50cm	60gsm	500 sheets
	WHCREPE060	white	60 x 60cm	60gsm	500 sheets
	WHCREPE075	white	75 x 75cm	60gsm	250 sheets
	WHCREPE090	white	90 x 90cm	60gsm	250 sheets
	WHCREPE100	white	100 x 100cm	60gsm	250 sheets
	WHCREPE1014	white	100 x 140cm	60gsm	125 sheets
	WHCREPE120	white	120 x 120cm	60gsm	125 sheets
	WHCREPE1315	white	130 x 150cm	60gsm	125 sheets

A




B

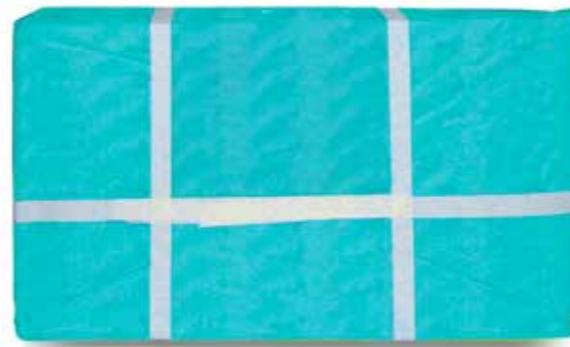

C



D

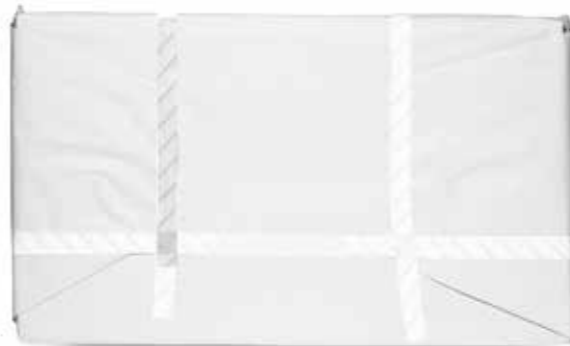

## Sterilisation Packaging – Liners



### A MILLENNIUM BLUE LINER

Lightweight and easy to use liners produced with a matte blue finish to minimise glare. 100% viscose-based product. 40gsm with similar absorbency to 110gsm tray liners. Compatible with steam sterilisation. Single use.

Product code	Size	Qty
MILB002	310 x 410mm	1000 sheets
MILB003	450 x 550mm	1000 sheets
MILB004	240 x 250mm	2400 sheets
MILB005	250 x 500mm	1200 sheets
MILB006	280 x 500mm	1200 sheets
MILB007	250 x 320mm	2000 sheets
MILB008	320 x 450mm	1200 sheets
MILB009	350 x 500mm	1200 sheets
MILB011	320 x 460mm	1000 sheets
MILB012	180 x 280mm	2000 sheets
MILB016	450 x 610mm	1000 sheets
MILB017	300 x 600mm	1200 sheets
MILB019	305 x 305mm	2000 sheets
MILB020	500 x 500mm	1000 sheets
MILB021	130 x 270mm	2000 sheets
MILB022	305 x 610mm	1000 sheets
MILB025	250 x 250mm	2400 sheets
MILB027	250 x 270mm	1000 sheets
MILB029	400 x 500mm	1000 sheets
MILB056	500 x 600mm	1000 sheets
MILB120X250	120 x 250mm	4800 sheets
MILB130X230	130 x 230mm	2000 sheets
MILB170X270	170 x 270mm	1000 sheets
MILB250X300	250 x 300mm	2000 sheets
MILB300X400	300 x 400mm	1000 sheets
MILB305X450	305 x 450mm	1200 sheets
MILB350X450	350 x 450mm	1200 sheets
MILB375X600	375 x 600mm	1000 sheets
MILB500X900	500 x 900mm	400 sheets



## Sterilisation Packaging – Liners



### A TRAY LINERS

Our Tray Liners have high absorbency and rapid dispersion to ensure condensation is removed during autoclaving, encouraging dry packs and improving the efficacy of the sterilisation process. To be used above or beneath trays to protect against basket. Thickness, 80gsm. Compatible with steam sterilisation. Single use.

Product code	Size	Qty
TLC001	305 x 15m	30 rolls/c
TLC002	457 x 15m	20 rolls/c
TLC004	200 x 305mm	2000 sheets
TLC005	250 x 250mm	2400 sheets
TLC006	250 x 500mm	1200 sheets
TLC007	305 x 410mm	1000 sheets
TLC008	457 x 550mm	800 sheets
TLC009	305 x 305mm	2000 sheets
TLC010	305 x 457mm	1000 sheets
TLC011	305 x 610mm	1000 sheets
TLC012	457 x 457mm	1000 sheets
TLC013	457 x 610mm	600 sheets
TLC014	250 x 490mm	1200 sheets
TLC017	400 x 457mm	1000 sheets
TLC018	305 x 540mm	1000 sheets
TLC019	430 x 457mm	1000 sheets
TLC020	350 x 457mm	1000 sheets
TLC022	360 x 540mm	1000 sheets
TLC023	356 x 457mm	1000 sheets
TLC024	457 x 692mm	600 sheets
TLC027	350 x 500mm	1000 sheets
TLC029	250 x 400mm	1200 sheets
TLC030	350 x 550mm	1000 sheets
TLC031	270 x 350mm	1000 sheets
TLC034	330 x 450mm	1000 sheets
TLC036	305 x 500mm	1000 sheets
TLC045	400 x 500mm	1000 sheets
TLC050	250 x 320mm	1000 sheets
TLC051	250 x 305mm	2000 sheets
TLC052	305 x 520mm	1000 sheets
TLC055	250 x 550mm	1200 sheets
TLC058	280 x 520mm	1000 sheets
TLC060	280 x 405mm	1000 sheets



A

K

Sterilisation Packaging





### TRAY LINERS CONTINUED

Our Tray Liners have high absorbency and rapid dispersion to ensure condensation is removed during autoclaving, encouraging dry packs and improving the efficacy of the sterilisation process. To be used above or beneath trays to protect against basket. Thickness, 80gsm. Compatible with steam sterilisation. Single use.

Product code	Size	Qty
TLC066	250 x 440mm	1000 sheets
TLC069	140 x 230mm	2000 sheets
TLC071	240 x 457mm	1000 sheets
TLC072	178 x 279mm	2400 sheets
TLC078	305 x 660mm	1000 sheets
TLC079	457 x 480mm	1000 sheets
TLC081	305 x 355mm	1000 sheets
TLC085	250 x 480mm	1000 sheets
TLC086	180 x 280mm	1000 sheets
TLC087	380 x 380mm	1000 sheets
TLC095	290 x 560mm	1000 sheets
TLC096	250 x 530mm	1200 sheets
TLC098	275 x 290mm	2000 sheets
TLC100	305 x 350mm	1000 sheets
TLC102	160 x 260mm	2400 sheets
TLC103	290 x 440mm	1000 sheets
TLC104	240 x 280mm	2400 sheets
TLC107	130 x 175mm	2000 sheets
TLC108	350 x 610mm	1000 sheets
TLC109	130 x 230mm	2000 sheets
TLC110	150 x 250mm	2000 sheets
TLC114	250 x 415mm	1000 sheets
TLC115	315 x 525mm	1000 sheets
TLC116	457 x 500mm	1000 sheets
TLC120	260 x 555mm	1200 sheets
TLC122	457 x 630mm	400 sheets
TLC220	220 x 305mm	2000 sheets
TLC3044	305 x 440mm	1000 sheets
TLC340	250 x 340mm	2000 sheets
TLC3535	350 x 350mm	1000 sheets
TLC3548	350 x 480mm	1000 sheets
TLC5060	500 x 600mm	600 sheets

