

INTRODUCING THE **NEW**
ELDER CARE ELECTRIC BED

WITH SAFETY FULL SIDE RAILS

FOR NURSING HOMES AND COMMUNITY HOSPITALS



DESIGNED TO BOOST

PRODUCTIVITY

STREAMLINE PATIENT HANDLING

IMPROVE **PATIENT SAFETY**

ELDER CARE ELECTRIC BED

FOR NURSING HOMES AND COMMUNITY HOSPITALS



POSITIONING AT YOUR FINGERTIPS

Handset controller with easy access for height adjustment and elevation



backrest



footrest



back & footrest



height elevation

ROBUST HEAD & FOOT BOARD

Removable head & foot board made from light weight ABS material with high pressure lamination panel colour of your choice. Option to secure the head or foot board to bed frame



IV POLE

Removable stainless steel fixed height IV pole with double hook. Two IV sockets at headend of bed locks IV pole securely.



EASY ACCESS CENTRAL BRAKE

Central lock pedal with 3 modes

- ✓ Total Lock
- ✓ Directional Steer
- ✓ Swivel

4-PANEL MATTRESS PLATFORM

Ventilated platform with three actuators for adjusting backrest (75°), footrest (30°) and height (340mm). Long lasting anti-rust framework treated with zinc-phosphate prior to epoxy powder coating.



ELDER CARE ELECTRIC BED

FOR NURSING HOMES AND COMMUNITY HOSPITALS

Framework Material	: Mild Steel
Finishing Treatment	: Zinc-Phosphate with Epoxy Powder Coating
Overall Dimension ⁽¹⁾	: L2040 x W955 x H360 -700 mm
Mattress Dimension ⁽²⁾	: L1860 x W870 x H100 mm
Sleeping Area ⁽³⁾	: L1910 x W860 mm
Height Adjustment	: 360 to 700 mm
Castor Wheel	: Ø125mm swivel castor with central lock
Nett Weight	: 92 kg
Safe Working Load	: 220 kg

Note

1. Overall length and width customisable
2. Mattress dimension customisable
3. Sleeping area may vary for customised overall dimension

SAFETY SIDE RAILS

Collapsible safety side rails with 3/4 side rails and 1/4 side rails on both sides of the bed to ensure patient safety. Bright colour side rails allow visual orientation to exit from bed easily. When side rails fully collapsed, it provides unobstructed access for lateral patient transfer.



MATTRESS

Medical grade water resistance "Medisoft" vinyl upholstered foam mattress for patient comfort.

WALL BUMPERS

Wall protection bumpers at four corners of bed

CASTOR WHEEL

Ø125mm swivel castors with central lock

URINE BAG HOOK

Convenient built-in double hook on both sides of bed to hang urine bag.




CENTRAL LOCK

Easy access central pedal for Locking all castor wheels (∨) Position for directional steer (–) Unlock all castor wheels to swivel (∧)



BED POSITION GUIDE

Backrest (75°), footrest (30°) and height elevation (340mm) with electric actuators

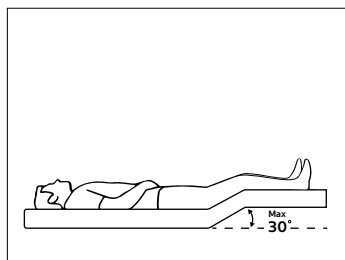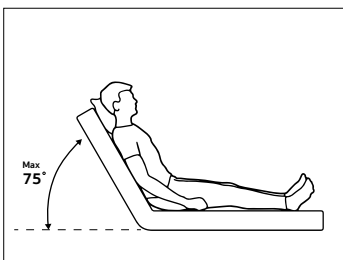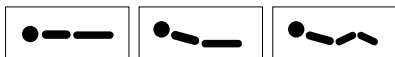

ORDERING INFORMATION

Brand : Meditron
Model : HB231UL-2CB-QCL-[Options]

COMPLIANCE

ISO13485 Medical Devices Quality Management Systems
IEC60601-1 Safety and Effectiveness of Medical Electrical Equipment





sales@neuindustries.com

 +65 9820 9495



www.neuindustries.com