


ELECTRIC 3-FUNCTION BED

Model V6k



V6k

FEATURES AND HIGHLIGHTS

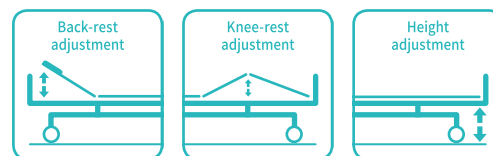

Technical parameters

▪ External size (LxW):	2085x995 mm
▪ Mattress platform:	1925x900 mm
▪ Safe working load:	200 kg
▪ Height adjustment:	420-730mm
▪ Back-rest adjustment:	0-75° (±5°)
▪ Knee-rest adjustment:	0-46° (±5°)

Optional

☐ Battery Back-up
☐ Bed Extension
☐ Nurse Controller
☐ Emergency Stop Button
☐ Auto-regression Function (Back Section)
☐ Night Light
☐ Manual CPR
☐ X-ray Translucent Mattress Platform (Back Section)
☐ Utility Shelf
☐ Trapeze Bracket/Pole

Technical configuration

☒ Linear motor	3 pcs
☒ Hand controller	1 pc
☒ 5" covered castors	4 pcs
☒ Foldable side rail	1 set
☒ PP bed ends with safe lock with Ø65 bumper wheels	1 set
☒ PP bed platform with integrated mattress retainer	1 set
☒ IV pole 4 hooks	1 set
☒ IV pole prevision	4 pcs
☒ Drainage hook	2 pcs

Optional color (MOQ100pcs- color changing)

RAL 210 710 30	RAL 5023	RAL 6019
RAL 3015	Pale oak	White oak
Red oak		

QUALITY CONTROL

SAIKANG medical, as a pioneer in the field of medical furniture, always choose branded suppliers strictly to ensure the safety and reliability of raw material. SAIKANG adopt robot welding, automatic epoxy coating, efficient assembly line, dedicated operator, rigorous quality control, only to guarantee SAIKANG's commitment on quality



Design

The unique germ-free design meets hospitals' strict infection control requirement. Seamless PP head and foot board and platform ensure the instant and sufficient sterilization.



Steels | welding

Baosteel from Fortune Global 500. Panasonic Robotic ensure 360° full smooth welding.



Plastic | Painting

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.
11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.



Electric Motor

American motor, with UL/ROHS/EN standard. Noise under 50dB, with lifetime 20,000 times movement.

Accessories



Platform

Integrated PP bed platform, adopting advanced manufacturing technology. 4-5 section platform, enables independent installation and easy maintenance.

Enlarge the seat area when back part is raised, providing more comforts for the patients.

Reduce the depression of the sponge at the seat panel.



Castors

5" covered casters. TPR tire no worn out after running 10 km.

Anti-winding hard shell cover, united forming without bolts.

Pass dynamic test: pass obstacles 40000 times, under the standard 200kg load of the whole bed.


Head & Foot Board

LG Chem fresh PP material, Plug in type.
Equipped with the bumpers wheel, diameter 65mm, preventing collision.




Mattress Retainer

Platform integrated mattress retainer at the both sides of the back and leg section, with a boarder of 33mm height.

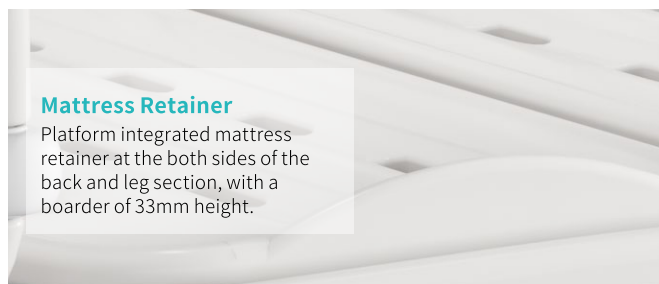


Guardrail

Aluminum alloy side rail, 6-staff, collapsed type with safe lock.

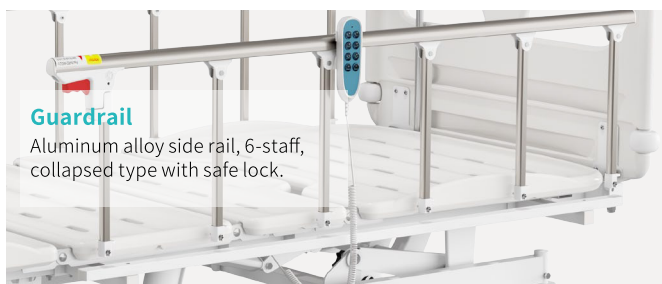

Wall bumper wheel

Horizontal corner bumpers at the outermost periphery of the bed to protect the bed and wall.

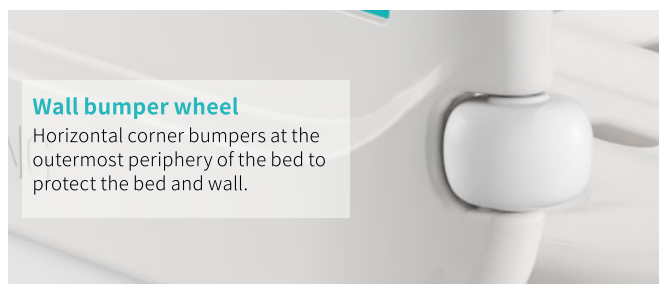

Drainage Hooks

Mobile plastic drainage hooks.




IV Pole

SS IV pole, bearing 15kg.



Back-rest adjustment

75°

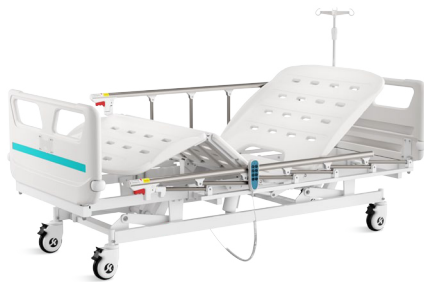

Knee-rest adjustment

46°



Height adjustment

420-730_{mm}



NEU Industries Pte Ltd

132 Joo Seng Rd, #04-01, Singapore 368358

Tel: +65 6817 8969 | Fax: +65 6636 9813

www.neuindustries.com | sales@neuindustries.com

