

FL250 Floorline Bed

Low Height, High Style



FL250 Floorline Bed



Performance Enhancing Accessories



BEA030050 - EasyMove Transport System

The EasyMove Transport System is an ideal transport solution designed to minimise manual handling. It allows effortless lifting and tilting of the bed to its transport position, while the large 100mm castors offer smooth and easy manoeuvring for one person. The EasyMove is also a perfect maintenance platform facilitating easy access to mechanical parts under the bed platform.



BEB046525 - Split folding Side Rails (Pair)

Left/Right pair of moulded rails provides manoeuvring assistance to a resident and/or a safety barrier. Can be positioned at the head end, foot end or a pair at both ends. Folds down to lowest position beside the bed platform.



BEB046520 - Folding Assist Bar

Provides manoeuvring assistance to a resident. Folds down when not in use, out of the way of transfers.



BEB046540 - Self Help Pole

The Self Help Pole provides manoeuvring and raising support to a resident. The Self Help Pole holder facilitates varied positioning and is securely attached behind the bed head.



BEB046530 - IV Pole

The height adjustable IV Pole can also be manoeuvred within its holder, which is discreetly attached behind the bed head.



BEB046590 - Battery Backup

Battery Backup is a recommended safety option to prolong bed operation by up to 15 lifts during a power outage.



BEB046500 - Width Extension Kit

Allows conversion of a Single FL250 into a King Single. Note - King Single FL250 Bed also sold separately.

Product Codes

- BEB046100 FL250 Floorline Bed - Single
- BEB046150 FL250 Floorline Bed - King Single

Bed End Sets - Includes Integrated Control Panel

- BEB046200 Bed End Set - Single - Tahoe (Sublime Teak)
- BEB046210 Bed End Set - Single - Vail (Milano Walnut)
- BEB046220 Bed End Set - King Single - Tahoe (Sublime Teak)
- BEB046230 Bed End Set - King Single - Vail (Milano Walnut)

Bed Side Panels - Also choose appropriate size attachment hardware.

- BEB046240 Bed Side Panels - Tahoe (Sublime Teak) - Hardware separate
- BEB046250 Bed Side Panels - Vail (Milano Walnut) - Hardware separate
- GSP391630 Bed Side Panels Attachment Hardware - Single Bed
- GSP391640 Bed Side Panels Attachment Hardware - King Single Bed

10 Year Warranty

Frame

3 Year Warranty

Electronics

1 Year Warranty

Handset & Cables



FL250 Bed Included Features

- 1 250kg SWL
- 2 110mm Low Height / 740mm Max Height
- 3 Single and King Single Sizes
- 4 Four Section Wipe Clean Mattress Platform
- 5 Powered Backrest from 0 - 60°
- 6 Powered Knee Break from 0 - 35°
- 7 Auto-Regression in Backrest and Knee Break for Shear-force Reduction
- 8 Trendelenburg / Reverse Trendelenburg
- 9 Handset and Attendant Control Panels - Backlit and Lockable
- 10 Underbed Lighting for Falls Prevention
- 11 Intelligent Wall Protection Buffer
- 12 Integrated Bed Extension to 100mm and 200mm
- 13 Mattress Retainers at Sides and Foot End (Extendable)
- 14 In-Line Central Locking Brakes
- 15 Stylish, Customisable and Bespoke Bed Ends and Sides
- 16 One Button Chair Positioning
- 17 One Button Optimum Ingress / Exit and CPR Height
- 18 Full Range of Performance Enhancing Accessories

FL250 Floorline Bed - Low Height, High Style

Technical Specifications


Standards Compliance

Aidacare FL250 Beds are designed and manufactured where required in compliance with International Medical Device Safety standard IEC60601-1, IEC60601-2-52. All electronics are designed for use in Acute Care environments and comply with all relevant Electrical Safety Standards.

Technical Specifications

Patient Surface - Single	2000mm x 900mm
Patient Surface - King Single	2000mm x 1060mm
Overall Length	2088mm to 2288mm
Length Extension	100mm / 200mm
Overall Width - Single	900mm
Overall Width - King Single	1060mm
Patient Surface Height Range	110mm to 740mm
Safe Working Load	250kg
Backrest Range	0 - 60 [degrees]
Knee Bend Range	0 - 35 [degrees]
Underbed Clearance	Variable. Minimum 30mm
Mattress Platform	4 Sections
Product Weight - Single	105kg
Product Weight - King Single	115kg
Castors	8 x Twin 75mm


Comprehensive user manual, troubleshooting guide, safety checklist, assembly instructions available.

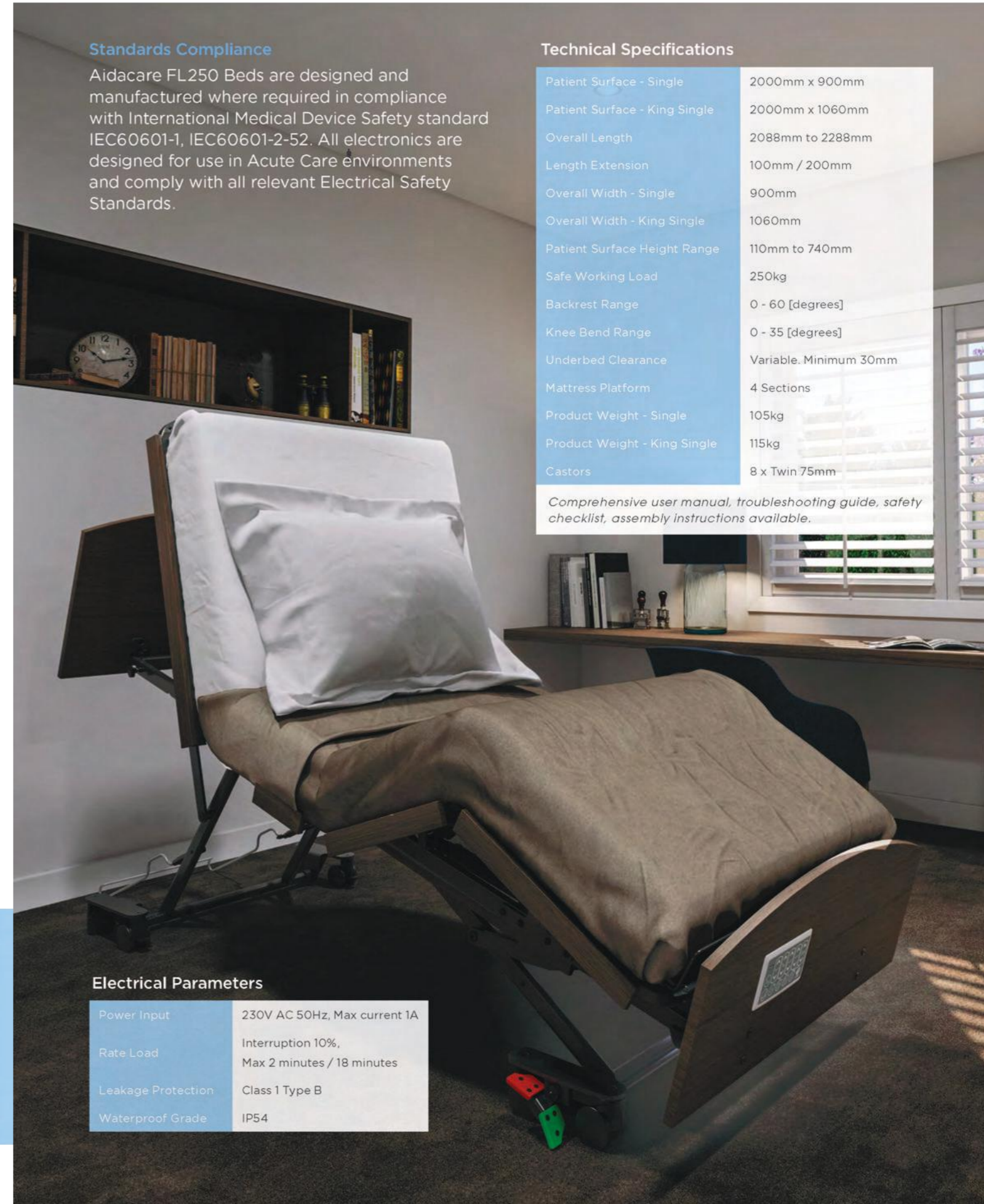

Electrical Parameters

Power Input	230V AC 50Hz, Max current 1A
Rate Load	Interruption 10%, Max 2 minutes / 18 minutes
Leakage Protection	Class 1 Type B
Waterproof Grade	IP54

Resident & Carer Safety

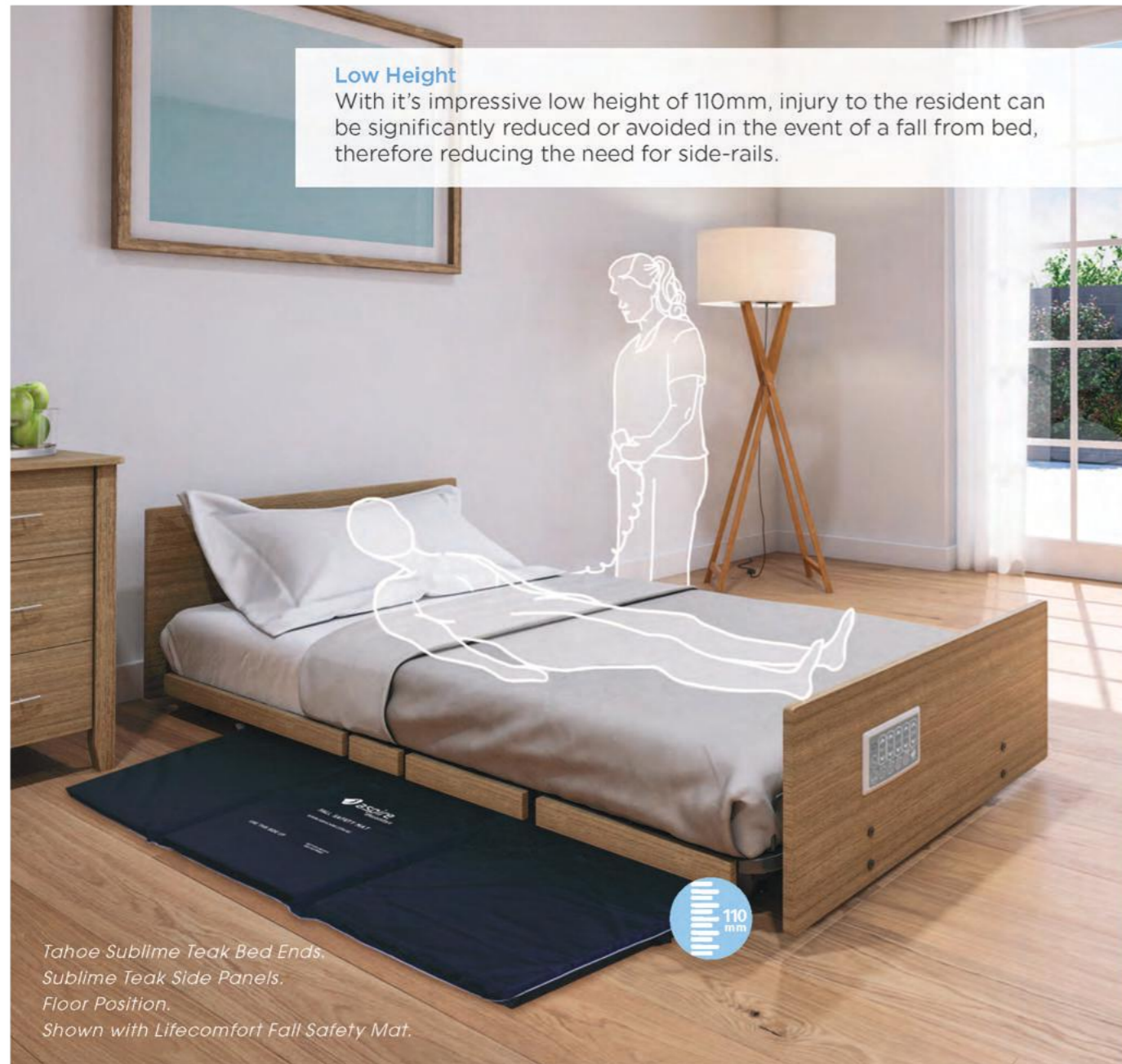


Attendant Control Panel

Allows the caregiver to lock out certain functions and also ensures bed operation is not just restricted to the handset.



Under Bed Clearance

While the FL250 achieves a low height of 110mm, full clearance under the bed is still maintained, minimising the risk of foot or cord impingement.

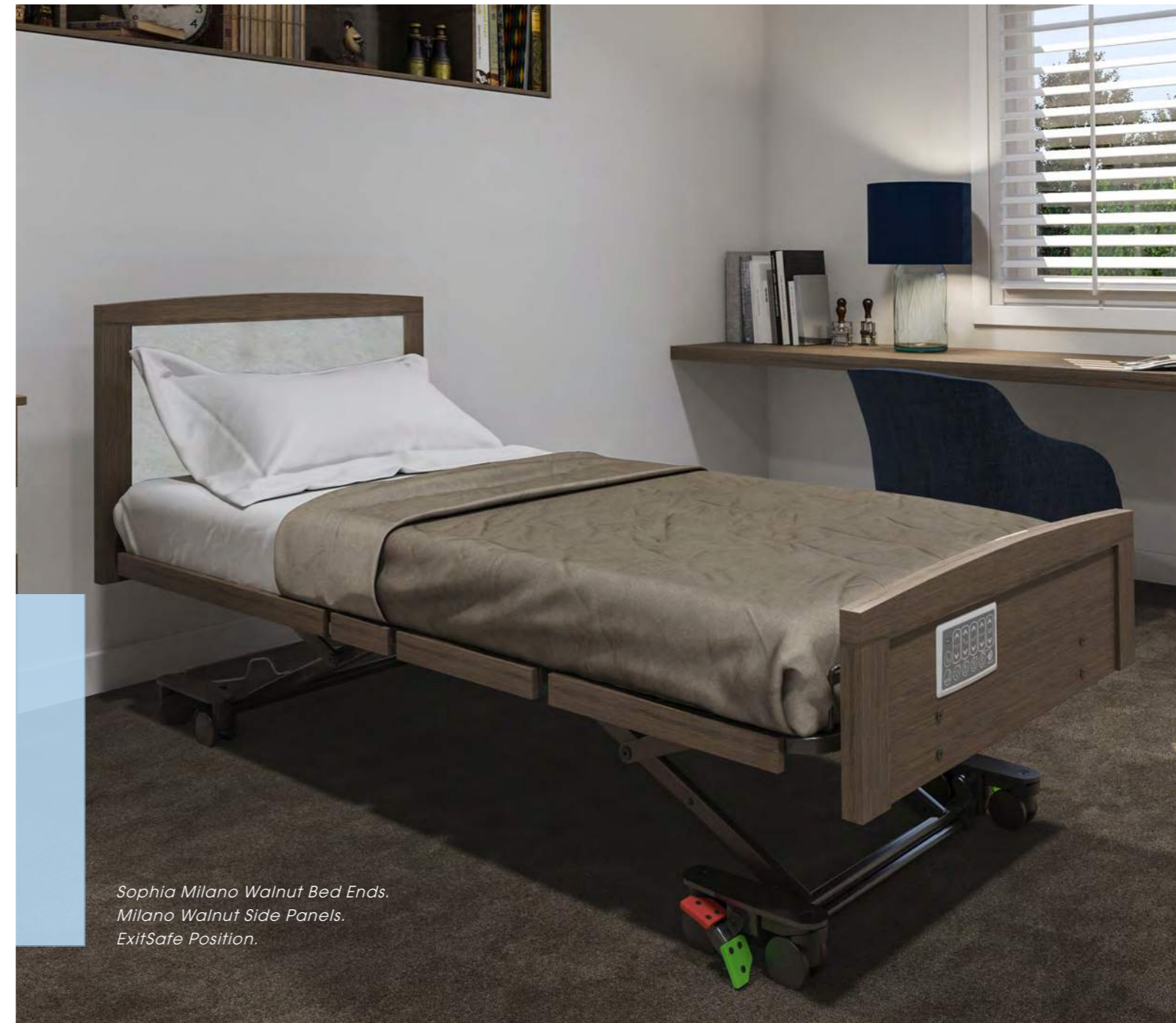

Customise to suit your interior design specifications with full split side rails or ABS end boards.



In Australian Aged Care Facilities, 1 in 2 residents fall each year. The bedside is the most common location for these falls.

1 Residents fall each year
in
2 INCREASED AGE
 HEALTH ISSUES
 COGNITIVE IMPAIRMENT

While falls can often be difficult to prevent, risk factors can be identified and strategies to reduce harm can be implemented. Sleeping low is a proven good strategy.

ExitSafe

By pushing the ExitSafe button the FL250 moves to a flat and ergonomic position at the optimum height for resident ingress / exit.



Lockable Functions



BrakeFirst

Brakes automatically lock when bed reaches it's lowest height. Brakes are located within the footprint of the bed, minimising any trip hazard for the caregiver and resident alike. Brakes are easily accessible at all heights.



SmartDescent

Automatically slows descent speed over the final 300mm of lowering to the floor level, allowing possible impingements to be identified and moved.



Underbed Lighting

Allows discreet visibility and safety at night.

Interior Design & Finishes


High Style

The FL250 Bed incorporates a sophisticated range of finishing options. From our contemporary stock selection and made to order styles through to completely bespoke, any interior design desires can be met.

Stock Bed Ends



Stock Side Panels

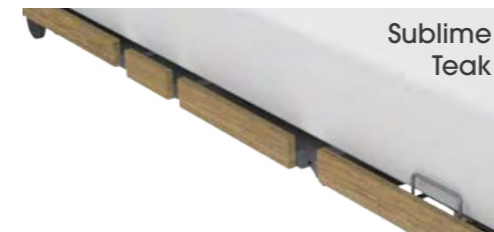

Made to Order Bed Ends



Tallo
 Timber or laminate surround and recessed timber veneer inner with a straight gap detail.



Eyre
 Timber or laminate surround and recessed timber veneer inner with a curved gap detail.



Como
 Timber or laminate surround and recessed timber veneer inner.



Sophia
 Timber surround with upholstered inner detail on headboard and recessed timber veneer detail on footboard.

Bed Side Panels are also available to match Made to Order Bed End colours and finishes. Materials may vary.

Bespoke Finishes

Don't see a Stock or Made to Order design that suits your needs? We can customise as required to incorporate your choice of material, colour, finish, size or shape and match any interior design specifications.

Clinical Considerations




Sleeping Low is Clinically Proven

International clinical studies have shown the benefits of sleeping low and eliminating side rails to resident safety, welfare and comfort, in particular by lowering fall rates and fall injury severity.

Introducing AutoRegression

AutoRegression helps to reduce sacral pressure, shear forces and abdominal compression during bed profiling.

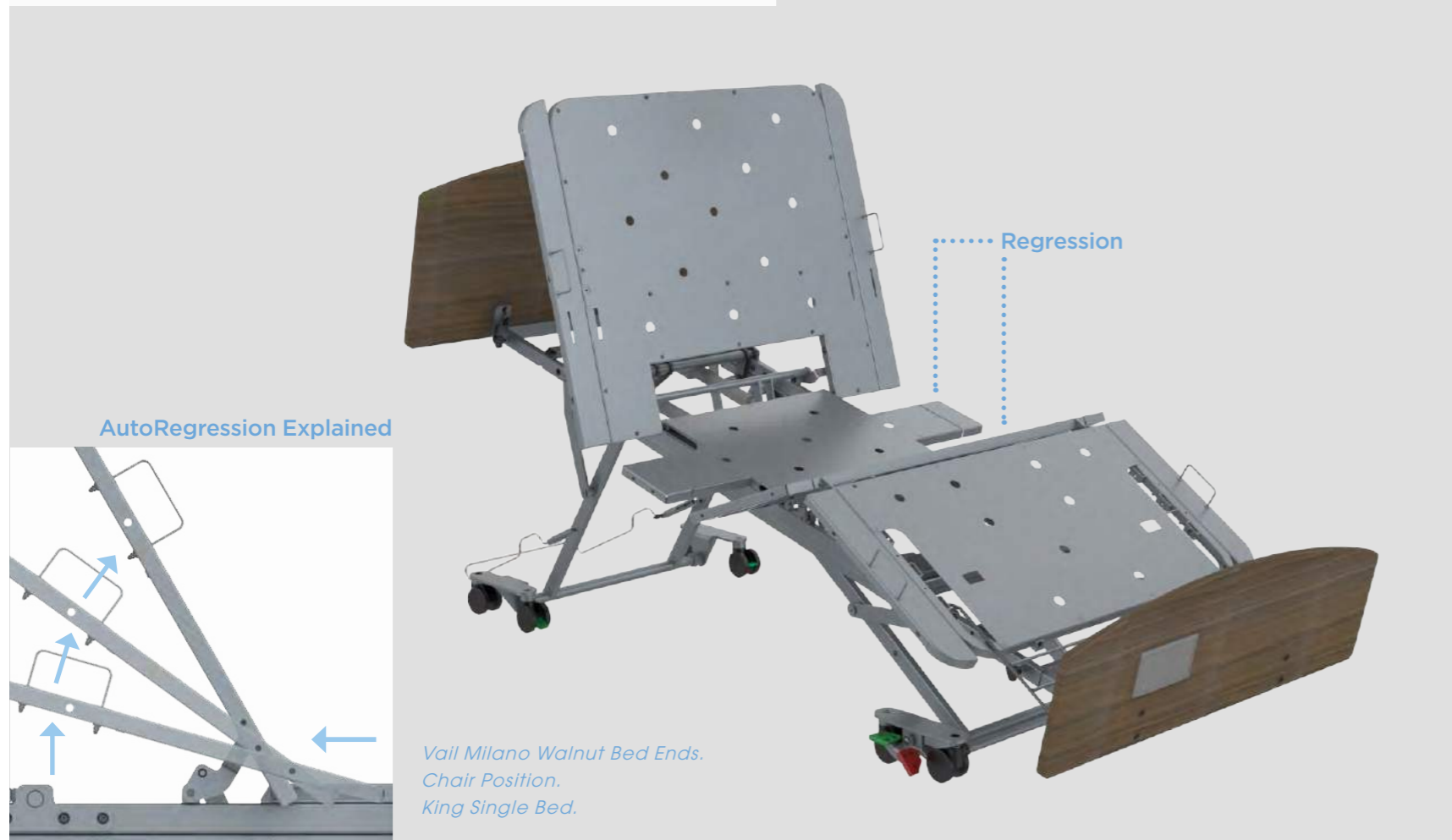

CPR Position

The low height and flat surface of the ExitSafe position also facilitates the administration of CPR in the event of an emergency.



Chair Positioning

For improved pressure management and respiration. Also a more comfortable dining and sitting position.




Full +/- Trendelenburg Functionality



Extendable Patient Surface

Integrated bed extension to 100mm or 200mm including foot end mattress retainer.

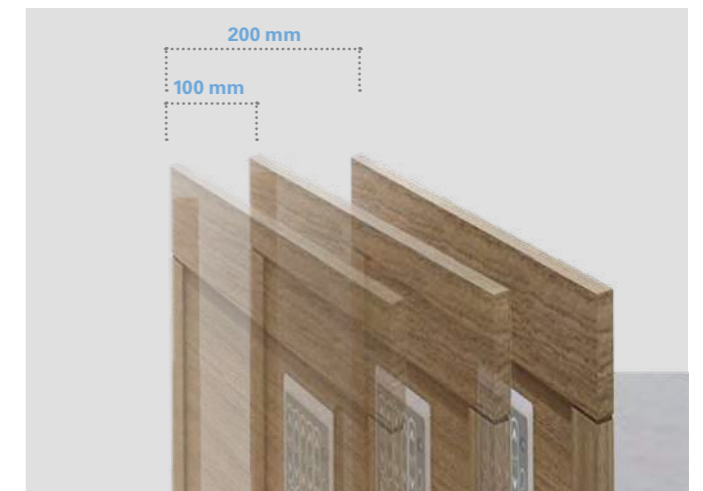

Integrating with Wireless Fall Prevention System



Accompanying Equipment & Furniture



ABB07.23 060418

Please speak to us about your made to order furniture requirements.
Visit our website to review our comprehensive range of equipment.



NEU Industries Pte Ltd
132 Joo Seng Road #04-01. Singapore 368358
Tel: +65 6817 8969 Fax: +65 6636 9813
Email: sales@neuindustries.com
www.neuindustries.com



Rev1(1222)