

FLOW 4 FC

Geriatric Recliner



Flow 4 FC

UPHOLSTERY

Self-extinguishing vinyl fabric (according to standard M2. EN1021.1-2) antibacterial and resistant to absorption and stains. Easy cleaning and maintenance. The chair upholstery is made of high-density polyurethane foam, designed for optimal comfort and ergonomic adaptation.

AVAILABLE COLORS



BLUE



BROWN



ARM REST ADJUSTMENT

Polyurethane armrests, height adjustable in five different positions by means of an integrated button.

CENTRAL
LOCK



FOUR Ø125MM
CASTOR

RETRACTABLE
FOOTREST TRAY

Flow 4 FC

TILTING

The combined tilting system of the seat keeps the patient suspended in “zero gravity” avoiding pressure on the lumbar area, pelvis and leg rest.



Adjustable, padded, removable headrest cushion for personalised comfort.




Structure made of steel tube with epoxy coating.



Independent gas spring adjustment for the backrest and leg rest with Bloc-o-Lift for smooth, secure positioning



Front access for lifts and dining tables.

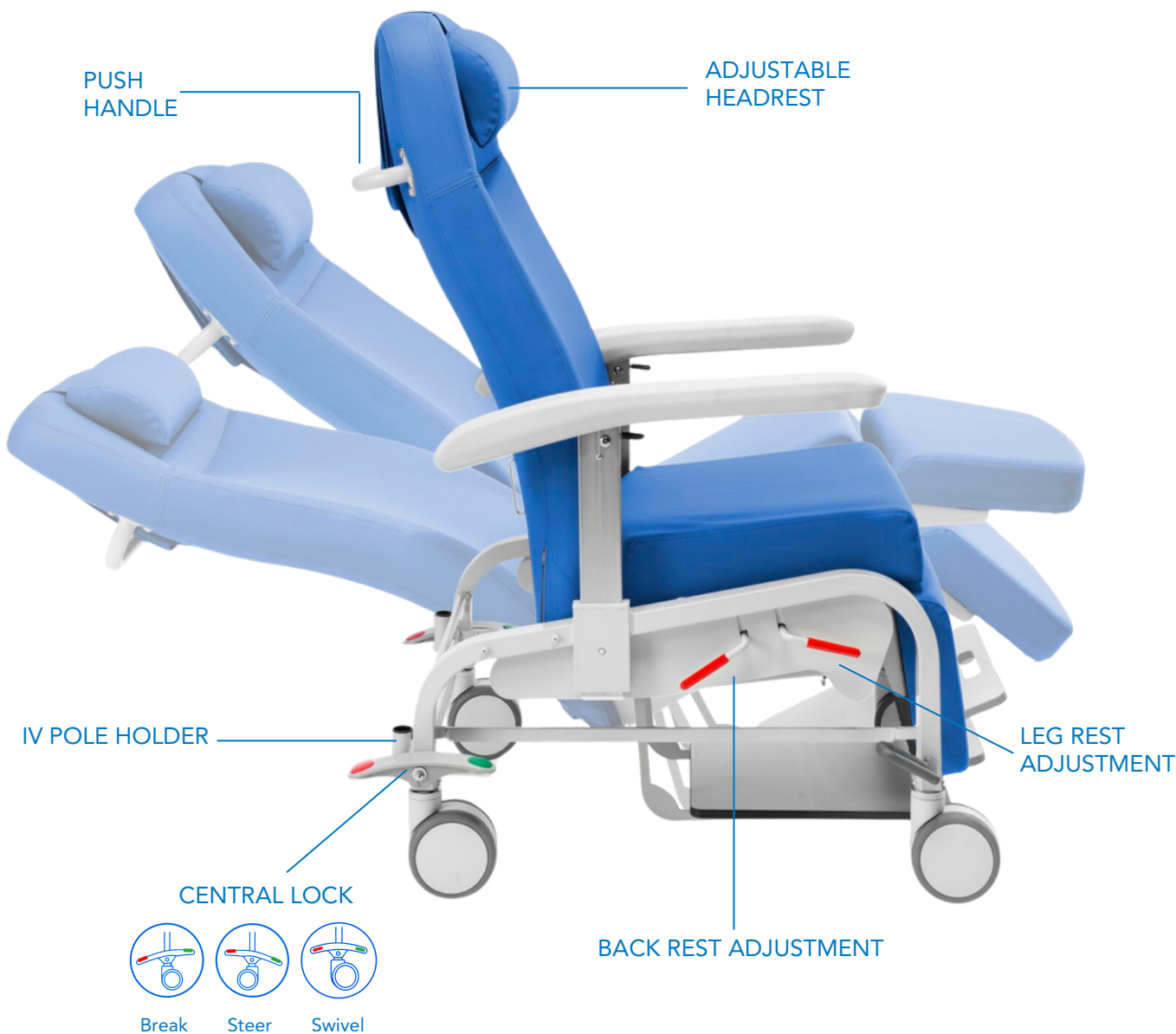

Flow 4 FC

TECHNICAL SPECIFICATIONS

Overall dimensions	H1225 × W720 × D860 mm
Extended length dimensions	1650 mm
Backrest dimensions	H840 × W510 mm
Seat dimensions	W520 × D520 × H460 mm
Legrest dimensions	W510 × D340 mm
Backrest angulation	85° × 168°
Extendable footrest	W290 - 350 × D220 mm
PUR armrest dimensions	W70 × D490 mm
Maximum load	180 kg
Approx. weight (no accessories)	45 kg



OPTIONAL ACCESSORIES



Seat belt



Meal Table



Oxygen Tank Holder



4-hook IV Pole



NEU Industries Pte Ltd

132 Joo Seng Rd, #04-01, Singapore 368358

Tel: +65 6817 8969 | Fax: +65 6636 9813

www.neuindustries.com | sales@neuindustries.com