



neu INDUSTRIES™
SAFER • CLEANER • BETTER

heel^{up}®

Pressure Relieving Heel Lift Cushion



LEVABO®
Simple solutions - Big impact

Content Page

Heel Up[®]

3

short

medium

long

max

Heel Up[®] **Air**

4

air medium

Heel Up[®] **Cut**

5

cut medium

Heel Up[®] **Fix**

6

fix short

fix medium

fix long

fix max

Heel Up[®] **Kids**

7

medium 15 cm

medium 20 cm

Accessories

8


Product data

Heel Up® is an environmentally friendly, inflatable, disposable positioning cushion. The product range has been developed in close cooperation with leading wound care nurses and ensures optimal comfort and hygiene in the prevention and treatment of pressure sores.





Danish developed
CE / MDR 2017
ISO 13485


short

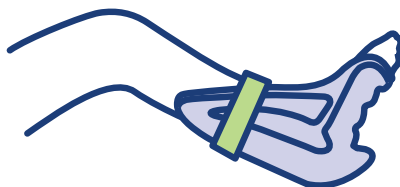


medium

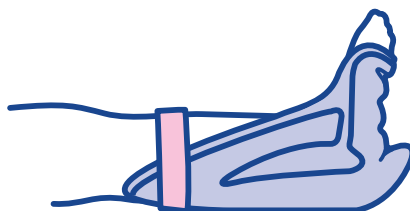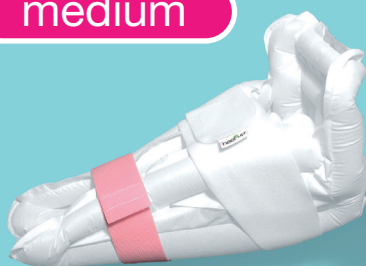

long

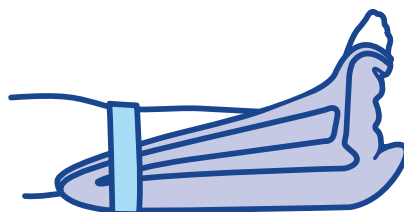

max

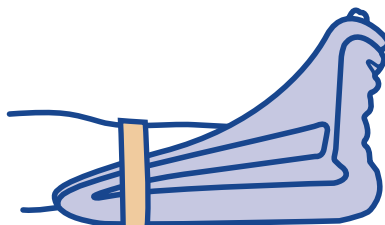

Product description

The heel lifter is designed to give complete heel relief. The heel lifter is shaped to fit the Achilles tendon, with height that supports the foot.

The ergonomic design of the heel lifter ensures stability without hindering the patient's movements.

Even with the slightest movement, the heel lifter will stimulate the microcirculation of the skin, which optimises circulation.

These properties make Heel Up® an optimal heel lifter when the bedridden patient needs a comfortable and pressure-relieving support surface.

Heel Up® is made of non-woven material, which makes the heel lifter soft, skin-friendly, heat and moisture absorbent and comfortable with optimal pressure relief.

Advantages

- Designed for the prevention and treatment of pressure ulcers in the heel region
- Single-use ensures high hygiene
- Easy to use. Heel Up® is inflated before use and adjusted to fit the individual user's anatomy and needs
- Good comfort. Soft, skin-friendly, heat and moisture-absorbent material
- The products eliminate the need for handling in a reuse depot and free up treatment resources
- Supplied in clean, sealed flat packaging including an inflation device
- Takes up minimal space in the storage area
- Tested and approved in accordance with ISO 10993-10:2010 skin irritation and sensitisation

Product name	Size	Packaging	Item No.
Heel Up Short	18 x 32 x 25 cm	12 pcs.	70001
Heel Up Medium	18 x 39 x 25 cm	12 pcs.	70007
Heel Up Long	18 x 48 x 25 cm	12 pcs.	70002
Heel Up Max	18 x 50 x 30 cm	12 pcs.	70018

Heel Up® pat. pending. Developed and produced in Denmark by Levabo ApS. Products are CE approved. Info is subject to change.

Heel Up® is an environmentally friendly, inflatable, disposable positioning cushion. The product range has been developed in close cooperation with leading wound care nurses and ensures optimal comfort and hygiene in the prevention and treatment of pressure sores.

ISO 13485
BUREAU VERITAS
Certification



Danish developed
CE / MDR 2017
ISO 13485

Product description

Heel Up® Air Medium is a heel lifter primarily intended for mobile patients. The heel lifter features an air cushion under the heel that helps prevent pressure injuries when the heel rests against a surface.


The ergonomic design of the heel lifter ensures stability without hindering the patient's movements.

Even with the slightest movement, the heel lifter will stimulate the microcirculation of the skin, which optimises circulation.

These properties make Heel Up® an optimal heel lifter when the bedridden patient needs a comfortable and pressure-relieving support surface.

Heel Up® is made of non-woven material, which makes the heel lifter soft, skin-friendly, heat and moisture absorbent and comfortable with optimal pressure relief.

Air Medium - The pressure relieving heel lifter

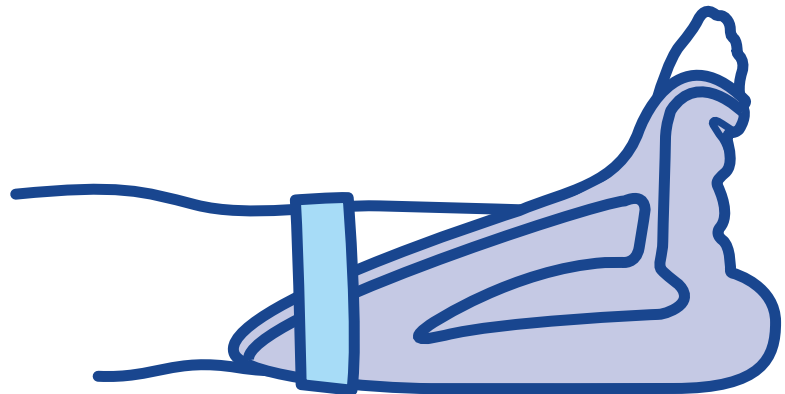

Advantages

- Designed for the prevention and treatment of pressure ulcers in the heel region
- Single-use ensures high hygiene
- Easy to use. Heel Up® is inflated before use and adjusted to fit the individual user's anatomy and needs
- Good comfort. Soft, skin-friendly, heat and moisture-absorbent material
- The products eliminate the need for handling in a reuse depot and free up treatment resources
- Supplied in clean, sealed flat packaging including an inflation device
- Takes up minimal space in the storage area
- Tested and approved in accordance with ISO 10993-10:2010 skin irritation and sensitisation



Product name	Size	Packaging	Item No.
Heel Up AIR Medium	18 x 39 x 25 cm	12 pcs.	70036

Product data

cut medium




Heel Up® is an environmentally friendly, inflatable, disposable positioning cushion. The product range has been developed in close cooperation with leading wound care nurses and ensures optimal comfort and hygiene in the prevention and treatment of pressure sores.



Danish developed
CE / MDR 2017
ISO 13485

Cut Medium - The pressure relieving heel lifter

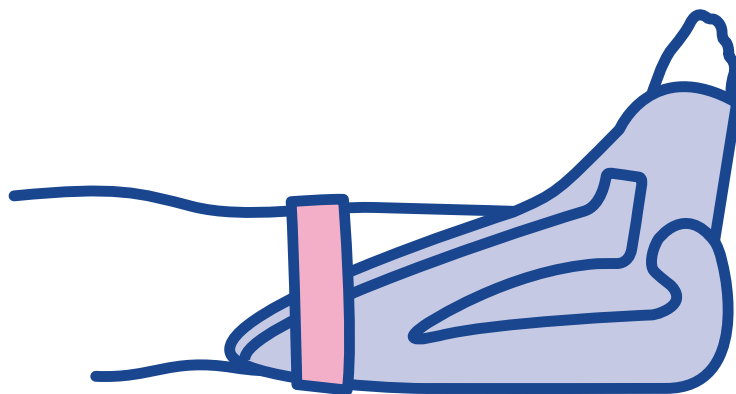

Product description

Heel Up® Cut Medium is a heel lifter primarily intended for patients with diabetes. The sail under foot is without air, allowing to make a cut-out for diabetic ulcers on the feet and toes.

The ergonomic design of the heel lifter ensures stability without hindering the patient's movements.

Even with the slightest movement, the heel lifter will stimulate the microcirculation of the skin, which optimises circulation.

These properties make Heel Up® an optimal heel lifter when the bedridden patient needs a comfortable and pressure-relieving support surface.

Heel Up® is made of non-woven material, which makes the heel lifter soft, skin-friendly, heat and moisture absorbent and comfortable with optimal pressure relief.

Advantages

- Designed for the prevention and treatment of pressure ulcers in the heel region
- Single-use ensures high hygiene
- Easy to use. Heel Up® is inflated before use and adjusted to fit the individual user's anatomy and needs
- Good comfort. Soft, skin-friendly, heat and moisture-absorbent material
- The products eliminate the need for handling in a reuse depot and free up treatment resources
- Supplied in clean, sealed flat packaging including an inflation device
- Takes up minimal space in the storage area
- Tested and approved in accordance with ISO 10993-10:2010 skin irritation and sensitisation

Product name	Size	Packaging	Item No.
Heel Up CUT Medium	18 x 39 x 25 cm	12 pcs.	70037

Heel Up® pat. pending. Developed and produced in Denmark by Levabo ApS. Products are CE approved. Info is subject to change.

Product data

Heel Up[®] is an environmentally friendly, inflatable, disposable positioning cushion. The product range has been developed in close cooperation with leading wound care nurses and ensures optimal comfort and hygiene in the prevention and treatment of pressure sores.



Danish developed
CE / MDR 2017
ISO 13485



Product description

Heel Up[®] Fix is a single-patient heel lifter designed to relieve pressure on the heel for mobile and seated users, particularly individuals in wheelchairs. The wide straps ensures that the heel lifter stays correctly in place, even if the user moves. Especially suitable for restless individuals.


The ergonomic design of the heel lifter ensures stability without hindering the patient's movements.

Even with the slightest movement, the heel lifter will stimulate the microcirculation of the skin, which optimises circulation.


These properties make Heel Up[®] an optimal heel lifter when the bedridden patient needs a comfortable and pressure-relieving support surface.

Heel Up[®] is made of non-woven material, which makes the heel lifter soft, skin-friendly, heat and moisture absorbent and comfortable with optimal pressure relief.


fix short

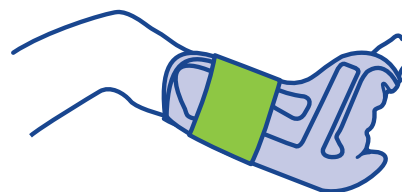


fix medium

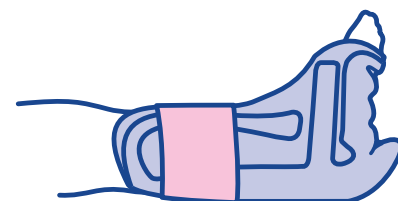

fix long

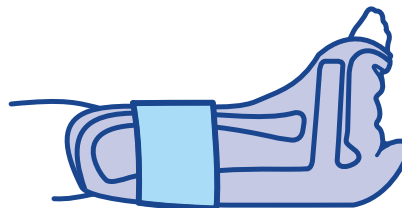

fix max

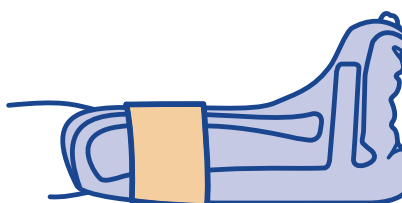

Advantages

- Designed for the prevention and treatment of pressure ulcers in the heel region
- Single-use ensures high hygiene
- Easy to use. Heel Up[®] is inflated before use and adjusted to fit the individual user's anatomy and needs
- Good comfort. Soft, skin-friendly, heat and moisture-absorbent material
- The products eliminate the need for handling in a reuse depot and free up treatment resources
- Supplied in clean, sealed flat packaging including an inflation device
- Takes up minimal space in the storage area
- Tested and approved in accordance with ISO 10993-10:2010 skin irritation and sensitisation

Product name	Size	Packaging	Item No.
Heel Up Fix Short	18 x 32 x 25 cm	10 pcs.	70146
Heel Up Fix Medium	18 x 39 x 25 cm	10 pcs.	70073
Heel Up Fix Long	18 x 48 x 25 cm	10 pcs.	70079
Heel Up Fix Max	18 x 50 x 30 cm	10 pcs.	70147

Heel Up[®] pat. pending. Developed and produced in Denmark by Levabo ApS. Products are CE approved. Info is subject to change.

Product data

Heel Up® is an environmentally friendly, inflatable, disposable positioning cushion. The product range has been developed in close cooperation with leading wound care nurses and ensures optimal comfort and hygiene in the prevention and treatment of pressure sores.




Danish developed
CE / MDR 2017
ISO 13485

Kids Medium - The pressure relieving heel lifter

medium 15 cm



15 cm

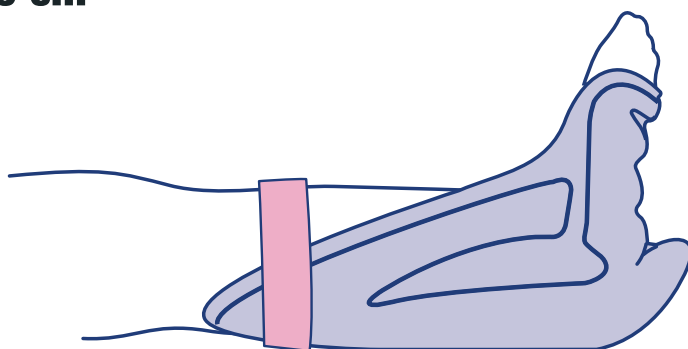


Lower Leg > 26 cm
Foot Length ≈ 15 cm
Appropriate for children aged 1½-4 years.

medium 20 cm



20 cm

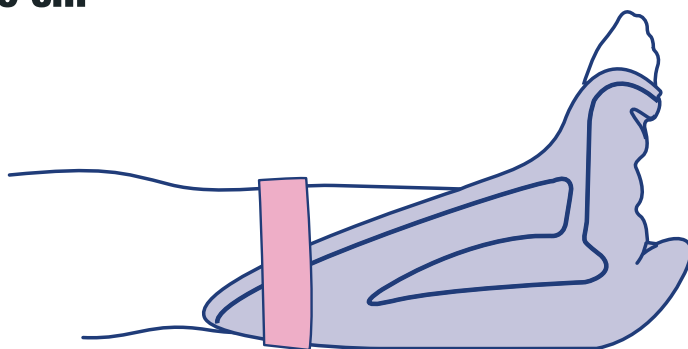

Lower Leg > 33 cm
Foot Length ≈ 20 cm
Appropriate for children aged 5-10 years.

Product description

Heel lifter for both prevention and treatment of pressure ulcers. Designed for bedridden children. Provides good support for the forefoot and minimizes the risk of foot drop.

The ergonomic design of the heel lifter ensures stability without hindering the patient's movements.

Even with the slightest movement, the heel lifter will stimulate the microcirculation of the skin, which optimises circulation.

These properties make Heel Up® an optimal heel lifter when the bedridden patient needs a comfortable and pressure-relieving support surface.

Heel Up® is made of non-woven material, which makes the heel lifter soft, skin-friendly, heat and moisture absorbent and comfortable with optimal pressure relief.

Advantages

- Designed for the prevention and treatment of pressure ulcers in the heel region
- Single-use ensures high hygiene
- Easy to use. Heel Up® is inflated before use and adjusted to fit the individual user's anatomy and needs
- Good comfort. Soft, skin-friendly, heat and moisture-absorbent material
- The products eliminate the need for handling in a reuse depot and free up treatment resources
- Supplied in clean, sealed flat packaging including an inflation device
- Takes up minimal space in the storage area
- Tested and approved in accordance with ISO 10993-10:2010 skin irritation and sensitisation

Product name	Size	Packaging	Item No.
Heel Up KIDS Medium 15 cm	11 x 26 x 16 cm	12 pcs.	70038
Heel Up KIDS Medium 20 cm	14 x 33 x 21 cm	12 pcs.	70039

Heel Up® pat. pending. Developed and produced in Denmark by Levabo ApS. Products are CE approved. Info is subject to change.

Accessories

Heel Up® is an environmentally friendly, inflatable, disposable positioning cushion. The product range has been developed in close cooperation with leading wound care nurses and ensures optimal comfort and hygiene in the prevention and treatment of pressure sores.



Danish developed
CE / MDR 2017
ISO 13485



240 V. AC Air Pump

Item Code: 76002

An electric pump that can be used for inflating the pressure relief cushions. Quiet operation and has an adjustable knob to control air pressure.

Air hose with valve included in packing.

Size: 24 x 13 x 11 cm

Packing: 1 pc

Power: 8 / 18 W

Pressure: 5-18 kPa

Output: 8 / 15 L/min



12 V. DC Battery Air Pump

Item Code: 76001

A portable battery-powered pump that can be used for inflating the pressure relief cushions. It runs quietly and is small and compact, easy to carry.

Size: 26 x 5 x 5 cm

Packing: 1 pc

Current: 375 mA (No Load)

700 mA (Max)

Pressure: 5-18 kPa

Output: 9 L/min



Manual Air Pump

Item Code: 76004

A small portable manual air pump that can be used for inflating the pressure relief cushions. Useful when no electricity or power is presently available. Its small and portable nature also means it is easy to carry around.

Size: 29 x 6 x 6 cm

Packing: 1 pc



Heel Up® Inlay

Item Code: 76015

Single-use non-woven felt inlay to be placed on a Heel Up® cushion to act as an additional protective layer for patients with high-sweating and to reduce risk of contamination to the cushion.

The inlay is placed on the cushion with the fabric side facing the patient's skin.

If necessary, the inlay can be cut into the required shape with a pair of scissors to ensure optimal air circulation and comfort around the wound area.

If the inlay is dirty or retains excessive moisture, it should be replaced.

Size: 35 x 50 cm

Packing: 10 pcs



NEU Industries Pte Ltd

132 Joo Seng Rd, #04-01, Singapore 368358

Tel: +65 6817 8969 | Fax: +65 6636 9813

www.neuindustries.com | sales@neuindustries.com