

TYPE 6 (Class 6) Emulating Indicator




Type 6 Emulating Indicators are part of an effective quality assurance program. It helps identify possible sterilization process failures that may result from improper packaging, improper loading or improper operation of the sterilizer. Emulating indicators change color when the specified parameter is provided. One emulator should be used in each sterilization package to ensure that the steam penetrates the packaging material and reaches the tools inside. These emulating indicators can be used for routine monitoring of steam sterilization cycles.


- * Water-based and non-toxic chemical indicator is used.
- * It is easy to read and interpret.
- * Compatible with ISO 11140-1.
- * Self-Adhesive versions can be used as a permanent record for traceability purposes.
- * Packed in 250 pieces.
- * Different designs are made for the company/person on request.




PRODUCT SIZE: 25 X 65 mm

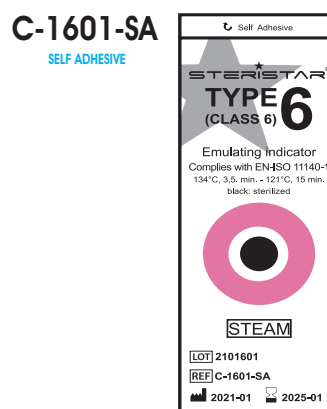

PRODUCT SIZE: 25 X 65 mm

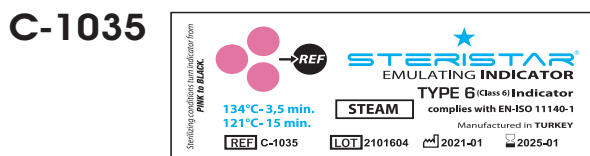

PRODUCT SIZE: 25 X 65 mm

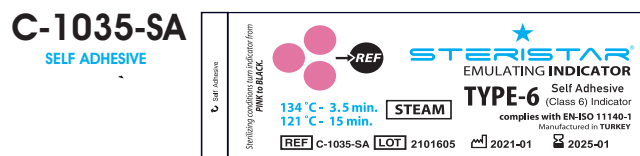

PRODUCT SIZE: 25 X 65 mm

REF NO	DESCRIPTION	PACKAGING
C-1035	Type 6 Emulating Indicator, Triple Indicators, 134 C-3,5 min.	80 Packs / CARTON
C-1035-SA	Type 6 Emulating Indicator, Triple Indicators, Self Adhesive, 134 C-3,5 min.	70 Packs / CARTON
C-1007	Type 6 Emulating Indicator, Triple Indicators, 134 C-7 min.	80 Packs / CARTON
C-1007-SA	Type 6 Emulating Indicator, Triple Indicators, Self Adhesive, 134 C-7 min.	70 Packs / CARTON
C-1601	Type 6 Emulating Indicator, Single Big Dot, 134 C-3,5 min.	80 Packs / CARTON
C-1601-SA	Type 6 Emulating Indicator, Single Big Dot, Self Adhesive, 134 C-3,5 min.	70 Packs / CARTON
C-1602	Type 6 Emulating Indicator, Single Big Dot, 134 C-7 min.	80 Packs / CARTON
C-1602-SA	Type 6 Emulating Indicator, Single Big Dot, Self Adhesive, 134 C-7 min.	70 Packs / CARTON